

A different approach to haemostasis.



Hemospray[®]
ENDOSCOPIC HEMOSTAT

For more information on Hemospray, visit
hemospray.cookmedical.com.

FOR
NONVARICEAL
**UPPER &
LOWER**
GI BLEEDS



Images courtesy of Dr John Morris,
Glasgow Royal Infirmary, Glasgow Scotland.





A different modality

What is Hemospray?

Hemospray is an inert mineral powder developed for endoscopic haemostasis that contains no human or animal proteins or botanicals and has no known allergens.


How does it work?

When Hemospray comes in contact with an actively bleeding site, the powder absorbs water, and then acts both cohesively and adhesively, forming a mechanical barrier over the bleeding site.

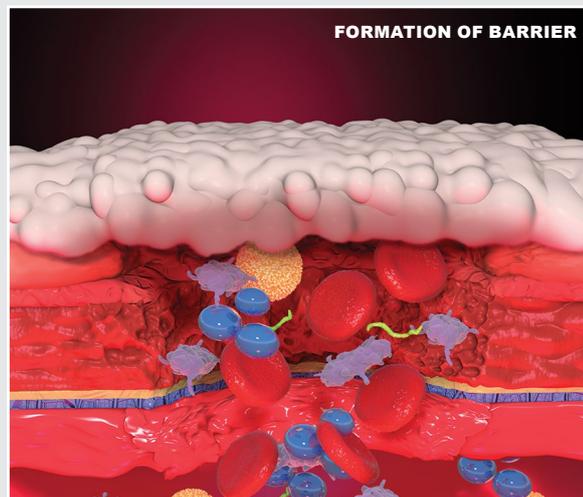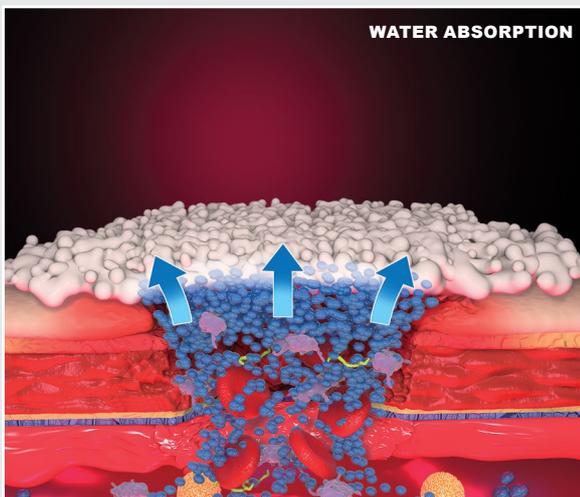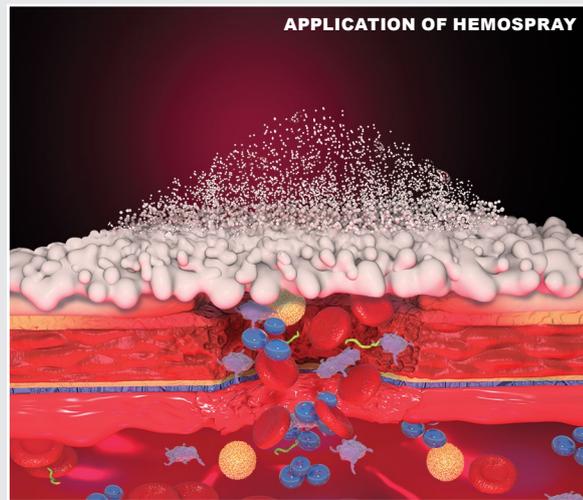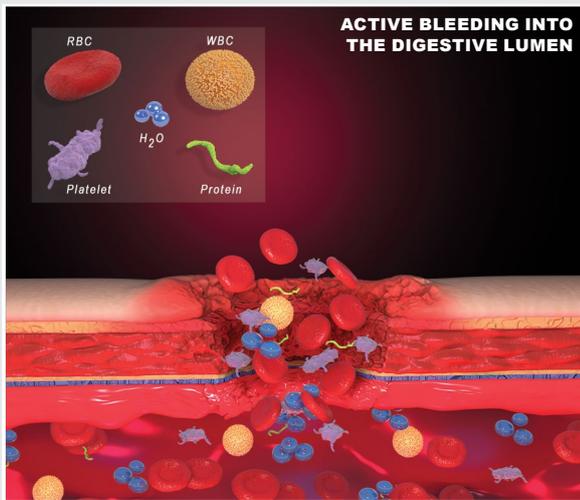

When to use Hemospray?

Hemospray is a different endoscopic modality that has demonstrated results in a wide range of nonvariceal gastrointestinal (GI) haemostasis procedures.

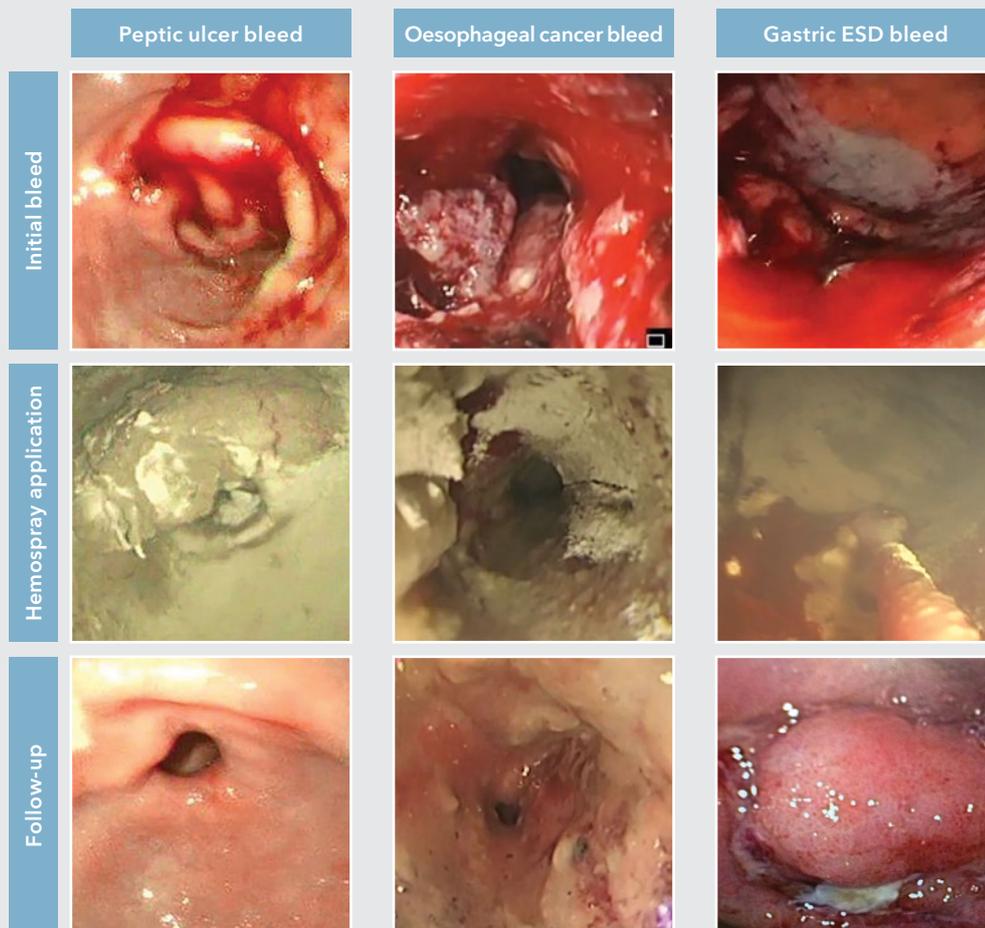

Hemospray has also demonstrated successful results with these additional nonvariceal bleed types:

Gastritis	Dieulafoy lesions
Mallory-Weiss tears	GAVE/watermelon stomach
Oesophagitis	Gastric angiodysplasia
GI neoplasms	Post endoscopic mucosal resection (EMR)

Peptic ulcer bleed images courtesy of Prof. Joseph Sung, Chinese University of Hong Kong, Hong Kong, China. Oesophageal cancer bleed image courtesy of Dr. John Morris, Royal Glasgow Infirmary, Glasgow, Scotland. Gastric ESD bleed images courtesy of Dr. David Serra, Hospital da Luz, Lisbon, Portugal.

A different approach to haemostasis.

Hemospray, unlike traditional therapies, is a nonthermal, nontraumatic, noncontact modality that doesn't require the precise targeting of other modalities.

- **Nonthermal:** The procedure causes no tissue changes, unlike some thermal modalities.
- **Nontraumatic:** Since no direct force is applied at the treatment site, the powder minimises the risk of tissue trauma that may result from other modalities.
- **Noncontact:** The aerosol delivery system eliminates the need for direct mechanical or contact forces with the treatment site by the delivery system.
- **Nonspecific targeting:** Powder is sprayed toward the source of the bleed, so it does not require the en face approach and the precise targeting of other modalities.

Simplifying haemostasis for you and your patient

Mechanical, injectable, and thermal therapies are used alone or in combination to manage GI bleeding. Hemospray is an easy to use, safe, and effective alternative to these traditional therapies. A systematic review of 20 studies¹ shows that Hemospray has a technical success rate of 97% when used to manage nonvariceal upper GI bleeding, with increasing success by publication year. Immediate haemostasis was achieved in 91% of patients. Immediate haemostasis, recurrent bleeding, and mortality were similar in studies comparing Hemospray to alternative methods. Complications were reported in <1% of patients.



technical success rate



achieved immediate haemostasis



patients with complications

1. Aziz M, Weissman S, Mehta TI, et al. Efficacy of Hemospray in non-variceal upper gastrointestinal bleeding: a systematic review with meta-analysis. *Ann Gastroenterol.* 2020;33(2):145-154.

Activate, open, and deploy.

In three steps you can deploy Hemospray to treat nonvariceal GI bleeds. Activate the CO₂ cartridge, open the valve, and spray toward the source of bleeding.



Hemospray reference articles

CLINICAL

Papaefthymiou A, Aslam N, Hussein M, et al. Hemospray® (hemostatic powder TC-325) as monotherapy for acute gastrointestinal bleeding: a multicenter prospective study. *Ann Gastroenterol*. 2024; 37(4): 418–426.

Pittayanon R, Khongka W, Linlaman S, et al. Hemostatic powder vs standard endoscopic treatment for gastrointestinal tumor bleeding: a multicenter randomized trial. *Gastroenterol*. 2023;165(3):762-772.e2.

Sung J, Moreea S, Dhaliwal H, et al. Use of topical mineral powder as monotherapy for treatment of active peptic ulcer bleeding. *Gastrointest Endosc*. 2022;96(1):28-35.

Lau J, Pittayanon R, Kwek A, et al. Comparison of a hemostatic powder and standard treatment in the control of active bleeding from upper nonvariceal lesions: a multicenter, noninferiority, randomized trial. *Ann Intern Med*. 2022;175(2):171-178.

Hussein M, Alzoubaidi D, Lopez M-F, et al. Hemostatic spray powder TC-325 in the primary endoscopic treatment of peptic ulcer-related bleeding: multicenter international registry. *Endoscopy*. 2021;53(1):36–43.

Hussein M, Alzoubaidi D, O'Donnell M, et al. Hemostatic powder TC-325 treatment of malignancy related upper gastrointestinal bleeds: international registry outcomes. *J Gastroentero Hepatol*. 2021. 2021;36(11):3027-3032.

Alzoubaidi D, Hussein M, Rusu R, et al. Outcomes from an international multicenter registry of patients with acute gastrointestinal bleeding undergoing endoscopic treatment with Hemospray. *Dig Endosc*. 2020;32(1):96–105.

Yen-I C, Wyse J, Lu Y, et al. TC-325 hemostatic powder versus current standard of care in managing malignant GI bleeding: a pilot randomized clinical trial. *Gastrointest Endosc*. 2020;91(2):321.e1–328.e1.

Hussein M, Alzoubaidi D, de la Serna A, et al. Outcomes of Hemospray therapy in the treatment of intraprocedural upper gastrointestinal bleeding post-endoscopic therapy. *United European Gastroenterol J*. 2020;8(10):1155–1162.

Hookey L, Barkun A, Sultanian R, et al. Successful hemostasis of active lower GI bleeding using a hemostatic powder as monotherapy, combination therapy, or rescue therapy. *Gastrointest Endosc*. 2019;89(4):865–871.

Mourad FH, Leong RW. Role of hemostatic powders in the management of lower gastrointestinal bleeding: a review. *J Gastroenterol Hepatol*. 2018;33(8):1445–1453.

Haddara S, Jacques J, Leclaire S, et al. A novel hemostatic powder for upper gastrointestinal bleeding: a multicenter study (the "GRAPHE" registry). *Endoscopy*. 2016;48(12): 1084–1095.

Chen YI, Barkun AN. Hemostatic powders in gastrointestinal bleeding: a systematic review. *Gastrointest Endosc Clin N Am*. 2015;25(3):535–552.

Chen YI, Barkun AN, Nolan S. Hemostatic powder TC-325 in the management of upper and lower gastrointestinal bleeding: a two-year experience at a single institution. *Endoscopy*. 2015;47(2):167–171.

Appleby VJ, Hutchinson JM, Beckett CJ, et al. Use of the haemostatic agent TC-325 in the treatment of bleeding secondary to endoscopic retrograde cholangiopancreatography sphincterotomy. *QJM*. 2015;108(1):79–80.

Masci E, Arena M, Morandi E, et al. Upper gastrointestinal active bleeding ulcers: review of literature on the results of endoscopic techniques and our experience with Hemospray. *Scand J Gastroenterol*. 2014;49(11):1290–1295.

Smith LA, Stanley AJ, Bergman JJ, et al. Hemospray for non-variceal upper gastrointestinal bleeding: results of the seal dataset (survey to evaluate the application of Hemospray in the luminal tract). *J Clin Gastroenterol*. 2014;48(10):e89–e92.

Sulz MC, Frei R, Meyenberger C, et al. Routine use of Hemospray for gastrointestinal bleeding: prospective two-center experience in Switzerland. *Endoscopy*. 2014;46(7):619–624.

Kratt T, Lange J, Königsrainer A, et al. Successful Hemospray treatment for recurrent diclofenac-induced severe diffuse lower gastrointestinal bleeding avoiding the need for colectomy. *Endoscopy*. 2014;46(suppl 1 UCTN):E173–E174.

Tarantino I, Barresi L, Granata A, et al. Hemospray for arterial hemorrhage following endoscopic ultrasound-guided pseudocyst drainage. *Endoscopy*. 2014;46(suppl 1 UCTN):E71.

Holster IL, Kuipers EJ, Tjwa ETTL. Hemospray in the treatment of upper gastrointestinal hemorrhage in patients on antithrombotic therapy. *Endoscopy*. 2013;45(1):63–66.

Chen YI, Barkun AN, Soulellis C, et al. Use of the endoscopically applied hemostatic powder TC-325 in cancer-related upper GI hemorrhage: preliminary experience (with video). *Gastrointest Endosc*. 2012;75(6):1278–1281.

NONCLINICAL

Shah ED, Law R. Valuing innovative endoscopic techniques: hemostatic powder for the treatment of GI tumor bleeding. *Gastrointest Endosc*. 2024;100(1):49–54.

Barkun AN, Adam V, Lu Y, et al. Using Hemospray improves the cost-effectiveness ratio in the management of upper gastrointestinal nonvariceal bleeding. *J Clin Gastroenterol*. 2018;52(1):36–44.

Beye B, Barret M, Alatawi A, et al. Topical hemostatic powder promotes reepithelialization and reduces scar formation after extensive esophageal mucosal resection. *Dis Esophagus*. 2016;29(6):520–527.

Product specifications


Used for haemostasis of nonvariceal gastrointestinal bleeding.

This device is supplied sterile and is disposable. It is intended for single use only.

Order Number	Reference Part Number	Catheter Fr	Catheter Length cm	Minimum Accessory Channel Diameter mm
G56572	HEMO-7	7.2	220	2.8
G21049	HEMO-10	9.6	220	3.7

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Distribution for details.

Device comes preassembled, and each package contains an extra catheter. No additional materials are required for usage. 3-year shelf life.

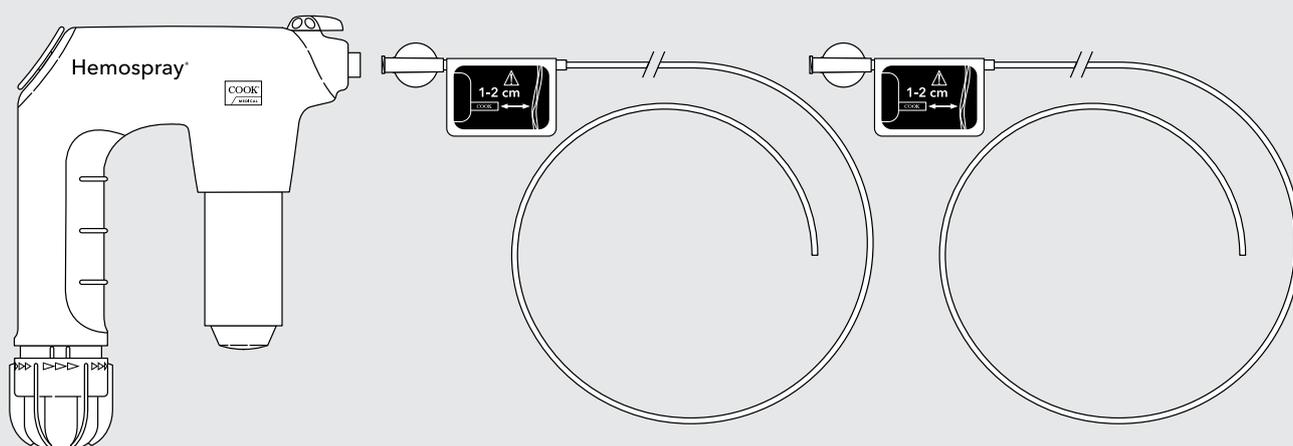

Refer to the current Instructions for Use (IFU) for detailed system use. Please see product risk information in the IFU at cookmedical.eu.

Customer Service

EU Website: cookmedical.eu
EDI: cookmedical.eu/edi
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001702, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Iceland: +354 800 7615, IS.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: orders@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: cau.custserv@cookmedical.com



AI, ESC, IR, OHNS, PI, RH, SUR-A4