


# A different approach to haemostasis.

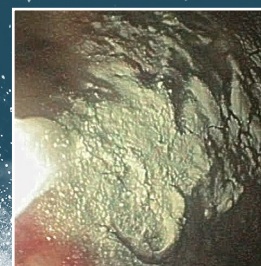

**Hemospray<sup>®</sup>**  
ENDOSCOPIC HEMOSTAT

For more information on Hemospray, visit  
[hemospray.cookmedical.com](http://hemospray.cookmedical.com).

FOR  
NONVARICEAL  
**UPPER &  
LOWER**  
GI BLEEDS



Images courtesy of Dr John Morris,  
Glasgow Royal Infirmary, Glasgow Scotland.



## Hemospray®

ENDOSCOPIC HEMOSTAT

This device is used for haemostasis of nonvariceal gastrointestinal bleeding.

This product is supplied sterile and is disposable. It is intended for single use only.

- **Nonthermal:** No tissue changes occur as sometimes experienced with thermal modalities.
- **Nontraumatic:** Since no force is applied at the treatment site, the powder minimises the risk of tissue trauma as experienced with other modalities.
- **Noncontact:** Aerosol delivery system eliminates the need for direct mechanical or contact forces with the treatment site by the delivery system.
- **Nonspecific targeting:** Powder is sprayed toward the source of the bleed, so it does not require the en face approach and the precise targeting of other modalities.

Order Number	Reference Part Number	Catheter Length cm	Catheter Fr	Minimum Accessory Channel Diameter mm
G56572	HEMO-7	220	7	2.8
G21049	HEMO-10	220	10	3.7

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Refer to the current Instructions for Use (IFU) for detailed system use.  
Please see product risk information in the IFU at [cookmedical.eu](http://cookmedical.eu).

### Customer Service

EU Website: [cookmedical.eu](http://cookmedical.eu)  
 EDI: [cookmedical.eu/edi](mailto:cookmedical.eu/edi)  
 Distributors: +353 61239240, [ssc.distributors@cookmedical.com](mailto:ssc.distributors@cookmedical.com)  
 Austria: +43 179567121, [oe.orders@cookmedical.com](mailto:oe.orders@cookmedical.com)  
 Belgium: +32 27001702, [be.orders@cookmedical.com](mailto:be.orders@cookmedical.com)  
 Denmark: +45 38487607, [da.orders@cookmedical.com](mailto:da.orders@cookmedical.com)  
 Finland: +358 972519996, [fi.orders@cookmedical.com](mailto:fi.orders@cookmedical.com)  
 France: +33 171230269, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Germany: +49 6950072804, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 Hungary: +36 17779199, [hu.orders@cookmedical.com](mailto:hu.orders@cookmedical.com)  
 Iceland: +354 800 7615, [is.orders@cookmedical.com](mailto:is.orders@cookmedical.com)  
 Ireland: +353 61239252, [ie.orders@cookmedical.com](mailto:ie.orders@cookmedical.com)  
 Italy: +39 0269682853, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Netherlands: +31 202013367, [nl.orders@cookmedical.com](mailto:nl.orders@cookmedical.com)  
 Norway: +47 23162968, [no.orders@cookmedical.com](mailto:no.orders@cookmedical.com)  
 Spain: +34 912702691, [es.orders@cookmedical.com](mailto:es.orders@cookmedical.com)  
 Sweden: +46 858769468, [se.orders@cookmedical.com](mailto:se.orders@cookmedical.com)  
 Switzerland - French: +41 448009609, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Switzerland - Italian: +41 448009609, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Switzerland - German: +41 448009609, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 United Kingdom: +44 2073654183, [uk.orders@cookmedical.com](mailto:uk.orders@cookmedical.com)

USA Website: [cookmedical.com](http://cookmedical.com)  
 EDI: [cookmedical.com/edi.do](mailto:cookmedical.com/edi.do)  
**Americas:**  
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335  
 E-mail: [orders@cookmedical.com](mailto:orders@cookmedical.com)  
**Australia:**  
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283  
 E-mail: [cau.custserv@cookmedical.com](mailto:cau.custserv@cookmedical.com)



AI-ESC-IR-OHNS-PI-RH-SUR-A4